


part number
CTP553

page 1 of 1

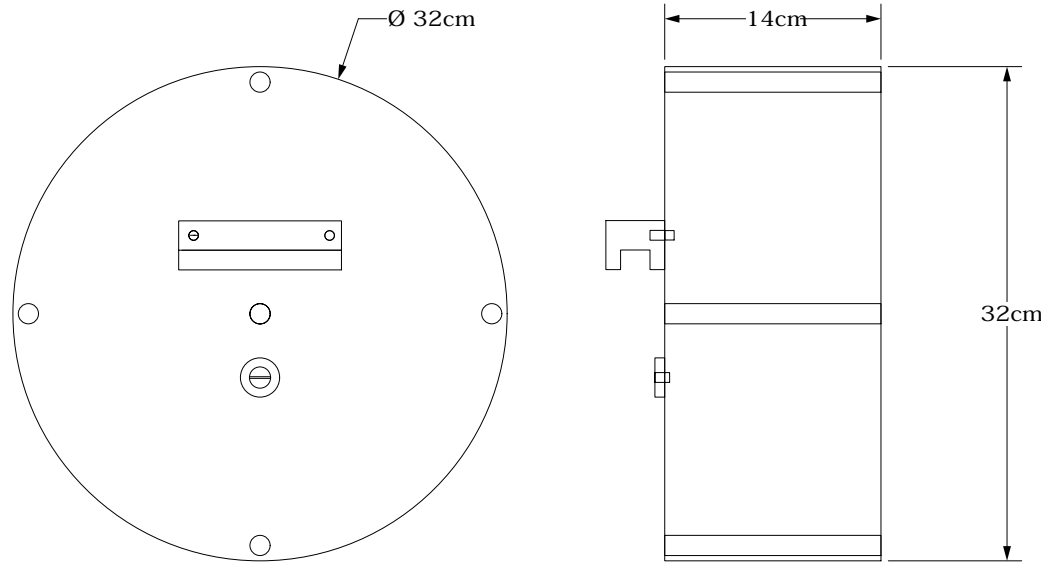

32cm Dose Phantom

The 32cm Dose Phantom is manufactured to comply with the US government's Center for Device and Regulatory Health (CDRH) dose phantom specifications. The specifications outline two different size phantoms, a 16cm diameter phantom for head dose evaluations and a 32cm diameter phantom for body dose evaluation. For specific data on the 16cm phantom review the CTP554 product information sheet.

The 32cm Dose Phantom is constructed out of acrylic. 13mm (1/2") holes are bored into the phantom providing cavities for any appropriately sized CT dose probes.

The phantom is supplied with 4 blank acrylic rods to be placed in the holes not used for the dose measurement. To ease mounting a case mount has been added to the phantom. This enables the phantom to be mounted on the standard Catphan® case.

The dose phantom is supplied with a wooden carrying case.



The Phantom Laboratory

phone
800-525-1190 or 518-692-1190
fax
518-692-3329

mail
P.O. Box 511
Salem, New York 12865
shipping
2727 State Route 29
Greenwich, New York 12834

email
info@phantomlab.com
web
www.phantomlab.com